



Class :3rd stage

Subject: Thermodynamics

م.م زهراء عبد الاله هادي

م.م ايلاف جاسم محان

م.م حسن علي شمخي

مهندس احمد عباس محمد



Ministry of Higher Education and Scientific Research

Al-Mustaqbal University College

Chemical engineering and petroleum industries (Thermodynamic Lab3)

Experiment No.6

Measurement of dryness Fraction by separating

Calorimeter

اعداد:

م.م زهراء عبد الاله هادي

م.م ايلاف جاسم محان

م.م حسن علي شمخي

مهندس احمد عباس محمد



Class :3rd stage

Subject: Thermodynamics

م.م زهراء عبد الاله هادي

م.م ايلاف جاسم محان

م.م حسن علي شمخي

مهندس احمد عباس محمد



Number of Experiment: Six

Name of Experiment: Measurement of dryness Fraction by separating calorimeter

Purpose of Experiment:

- 1- Measurement of dryness fraction of steam.
- 2- Understanding various methods of measurement of dryness fraction.
- 3- Arrangement of various components for set up of throttling and separating Calorimeter.

Equipment:

1. Boiler.
2. Separating calorimeter.
3. Throttling calorimeter.
4. Class container.
5. Pressure gauge.
6. Thermometer.
7. Condenser.

Theory:

Throttling and Separating Calorimeters are used to measure the dryness fraction of steam. If the steam is sufficiently dry to enable it to reach a supersaturated state by throttling process, then a solitary throttling calorimeter can serve the purpose. If on the other hand, the steam is very wet, then a separating calorimeter is used.

When the dryness fraction of the steam is somewhere in between, then a combined setup of Throttling and Separating Calorimeter is used

1. Separating calorimeter

The equality of wet steam is usually defined by its dryness fraction. When the dryness fraction, pressure and temperature of the steam are known, then the state of wet steam is fully defined. In a steam plant it's at times necessary to know the state



Class :3rd stage

Subject: Thermodynamics

م.م زهراء عبد الاله هادي

م.م ايلاف جاسم محان


م.م حسن علي شمخي

مهندس احمد عباس محمد



of steam. For wet steam, this entails finding the dryness fraction. When the steam is very wet, we make use of a separating calorimeter.

Construction of separating calorimeter is as shown in figure:

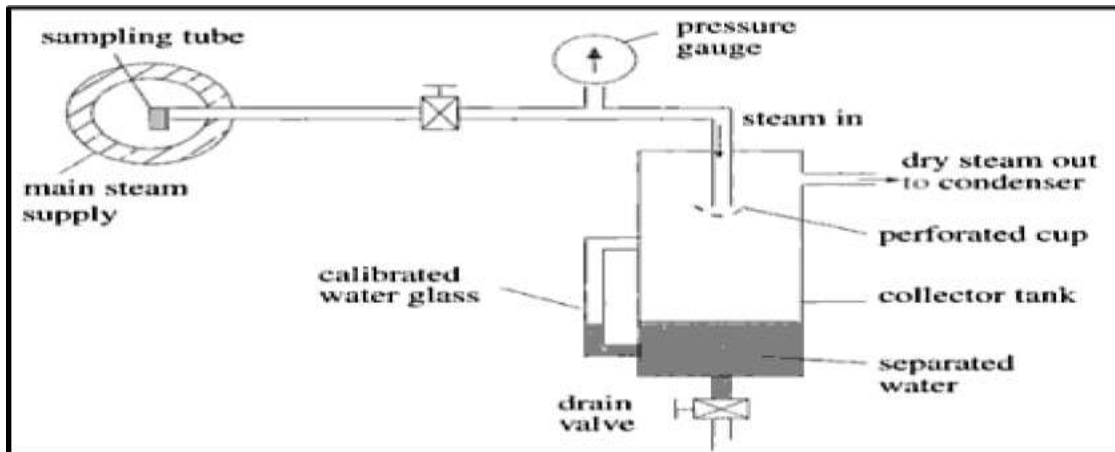

The steam is collected out of the main steam supply and enters the separator from the top. The steam is forced to make a sharp turn when it hits the perforated cup (or any mechanism that produce the same effect).this result in vortex motion in the stem and, and water separates out by the centrifugal action. The droplets then remain inside the separator and are collected at bottom, where level can be recorded from the water glass. The dry steam will pass out of the calorimeter into a small condenser for the collection of the condensate. However, not all the water droplets remain in the collector tank. Some water droplets pass through to the condenser, and hence this calorimeter only gives a close approximation of the dryness fraction of the steam.

From the results obtained from the two collectors, the dryness fraction may then be found from

2. Throttling calorimeter

If we have steam that is nearly dry, we make use of a throttling calorimeter as shown in Fig (2). This calorimeter is operated by first opening the stop valve fully



Class :3rd stage

Subject: Thermodynamics

م.م زهراء عبد الاله هادي

م.م ايلاف جاسم محان

م.م حسن علي شمخي

مهندس احمد عباس محمد




so that the steam is not partially throttled as it passes through the apparatus for a while to allow the pressure and temperature to stabilize. If the pressure is very close to atmospheric pressure, the saturation should be around 100 Co, it may be assumed that the steam is superheated.

When the condition has become steady, the gauge pressure before throttling is read from the pressure gauge. After throttling, the temperature and gauge pressure is read from the thermometer and manometer respectively. The barometric pressure is also recorded

From the equation: -

Dryness fraction at throttling calorimeter =

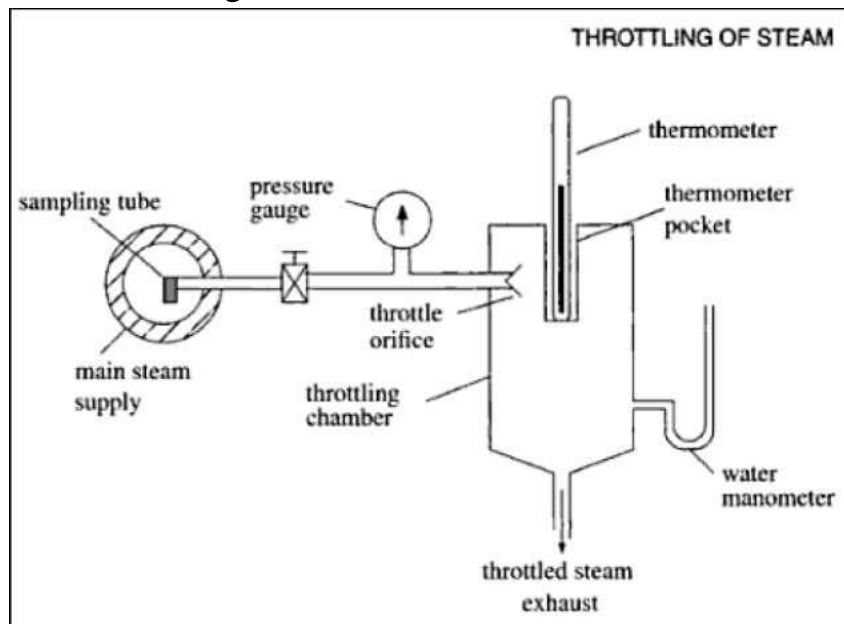

Procedure:

Before starting the experiment, it is necessary to remove the condensate water remaining from the previous experiment from the machine to make the calculations accurate. Use Figure 1 when performing the experiment steps as follows

1. Make sure that there is water inside the manometer scale at zero reading
2. Close the valve 1
3. Operate the steam boiler to the required pressure and make sure the steam



Class :3rd stage

Subject: Thermodynamics

م.م زهراء عبد الاله هادي

م.م ايلاف جاسم محان

م.م حسن علي شمخي

مهندس احمد عباس محمد



4. reaches the valve1
5. Open the valves 2-9
6. Close the terminals 10 and 11
7. Open valve 1 gradually to allow for steam entry and ensure that all condensation
8. residues are expelled from the previous test of LEDs 2 and 9
9. After a short while, the steam will flow into the first part of the oil. Open the
- 10.throttle valve 11 slowly and wait for a few minutes until the second throttle
- 11.reaches the stable temperature.
- 12.Close the valve 1 and wait for a few moments until the meter reading becomes zero
- 13.Shut the valves 2, 4, 5, 6, 7, 9, 11
- 14.Open the valves 3, 8, 10
- 15.Operate the temperature measuring device
- 16.Place a graduated baker with a volume greater than 2 liters below the 3 valve to
- 17.collect the condensate water. Place another vessel that is the same size as the
- 18.first one. Place it under valve 2 to collect the collected water from the separation
- 19.valve with the remaining 2, 3 closed
- 20.Then slowly open the throttle valve 11 and record the reading of the manometer
- 21.(P P_3) and try to control the opening of the throttle valve 11 to get stable
- 22.readings for several minutes and then record the reading of the pressure gauge
- 23.(P_2) And temperature (t_3).



Class :3rd stage

Subject: Thermodynamics

م.م زهراء عبد الاله هادي

م.م ايلاف جاسم محان

م.م حسن علي شمخي

مهندس احمد عباس محمد



24. Close the valve 1 and wait a little while to complete the water pool in the opening
25. price and then open the 2, 9 and measure the volume of water out of the V_1 and
26. the volume of the water out of the throttle (V_2)
27. Start calculating the dryness of wet steam



Discussion:

1. Mention the reasons for error in the results?
2. What is calorimeter?